

PLATE BAGS STIRRER

MOD. 723

COD. 29970018



DESCRIPTION

The Plate Bag Stirrer mod. 723 is equipped with a separate cell holder which can contain 18 bags cc. 400.

A constant speed of 5 revolutions per minute keeps the cellular material in continuous suspension before the transfusion.

The perforated basket allows air circulation so the temperature variations are avoided. The motor position guarantees continuous operation for several hours without overheating.

Outer casing in powdered painted steel acids resistant, basket made of stainless steel Aisi 304 2B glazed.

Main switch bright.

OPERATING FEATURES


TECHNICAL SPECIFICATIONS

- ♦ Work surface dimensions - 9 bins, WxDxH: 250 x 20 x 150 mm.
- ♦ Overall dimensions WxDxH: 520 x 320 x 360 mm.
- ♦ Weight: 10 Kg

ELECTRICAL SPECIFICATIONS

- ♦ Power supply: 230 V - 50 Hz.
- ♦ Electrical input: 45 W
- ♦ Overload fuses: 2 x 0.5 A (5x20) mm.
- ♦ Connection outlet: 10 A

PARTS LIST

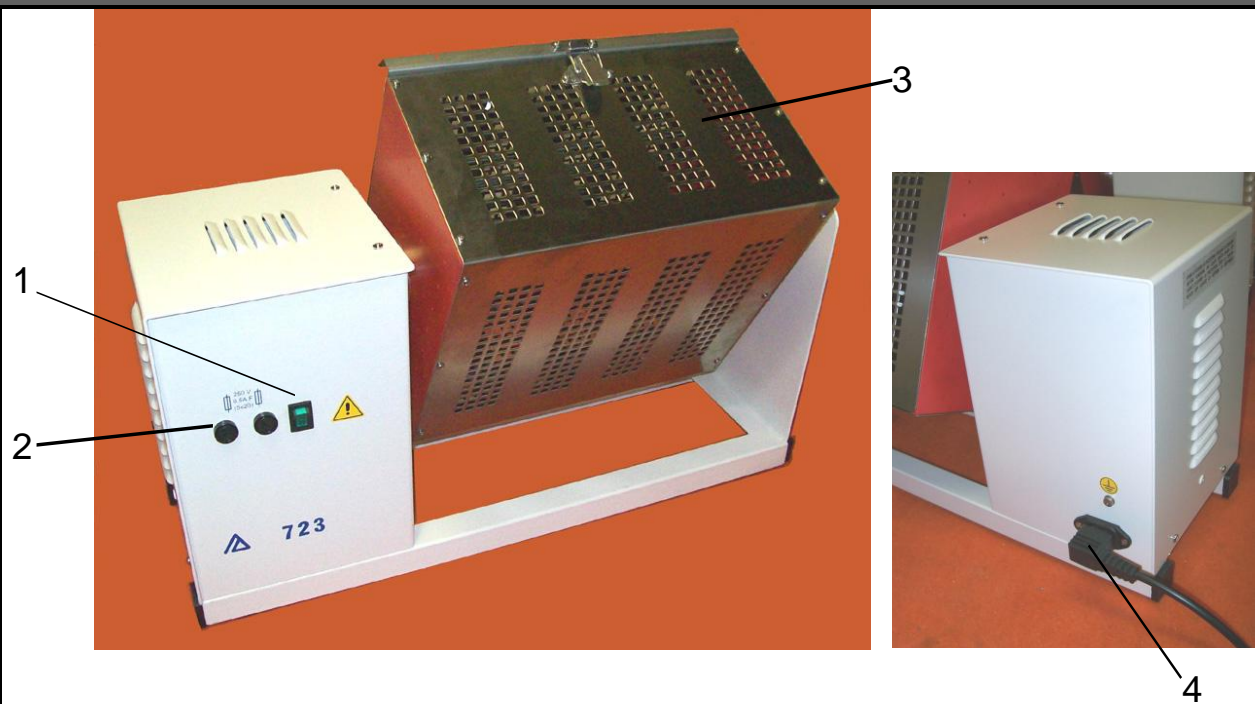

- 1 Main switch bright
- 2 Overload fuses
- 3 Cell holder basket
- 4 Power cable outlet