




VERTICAL SHAKING STIRRER

MOD. 709/V

COD. 29970017

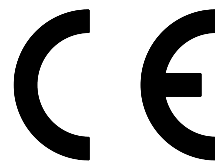

DESCRIPTION

The vertical shaking stirrer mod.709 / V has a working surface with a vertical movement of 30 mm, suitable for any type of agitation, and in particular in the extractions.

Electronic regulation of speed from 100 to 400 rpm (regulation 1/min, 10/min or 100/min).

Possibility of continuous or pre-set working from 1' to 999' minutes.

Possibility to memorise 4 disjoint programming channels of time and speed values.

Polycarbonate anti-scratch panel with visualisation of speed (RPM) and remaining time to finish cycle.

Outer casing in powdered painted steel acids resistant and work surface in stainless steel AISI 304 2B glazed.

Main switch bright.

The shaking vertical stirrer mod. 709/V needs a support for separating funnels, select from the following, or universal type basket, to be ordered separately as an "accessory":

- ◆ support for nr. 8 funnels of cc.100 or cc.250
- ◆ support for nr. 4 funnels of cc.500/1000 or nr. 2 funnels of cc. 2000 (do not apply more than two funnels of cc. 2000)
- ◆ universal type baskets in stainless steel AISI 304 2B glazed suitable to hold flasks, beakers etc.

OPERATING FEATURES


TECHNICAL SPECIFICATIONS

- ◆ Overall dimensions WxDxH: 630 x 455 x 665 mm.
- ◆ Working surface dimensions WxD: 550 x 330 mm.
- ◆ Weight: 33 Kg. with vertical support
- ◆ Max capacity of work surface: 5 Kg.

ELECTRICAL SPECIFICATIONS

- ◆ Power supply: 230 V - 50 Hz.
- ◆ Electrical input: 70 W
- ◆ Overload fuses: 2 x 2 AF (5 x 20) mm.
Connection outlet: 10 A

ACCESSORIES - SUPPORT FOR SEPARATORY FUNNELS

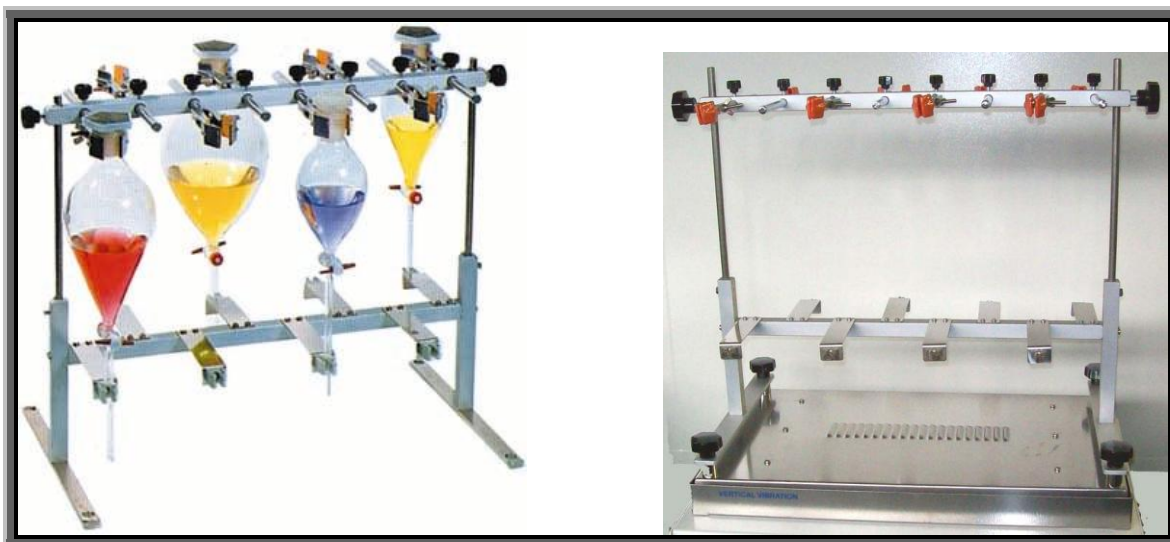

OPTION

Vertical stirrer mod. 709/V with universal basket suitable for holding and fixing, beakers and flasks of different sizes.

