




WELL D-ONE® SYSTEM

OYRON WELL D-ONE®

For Urinary tract infections (UTI)

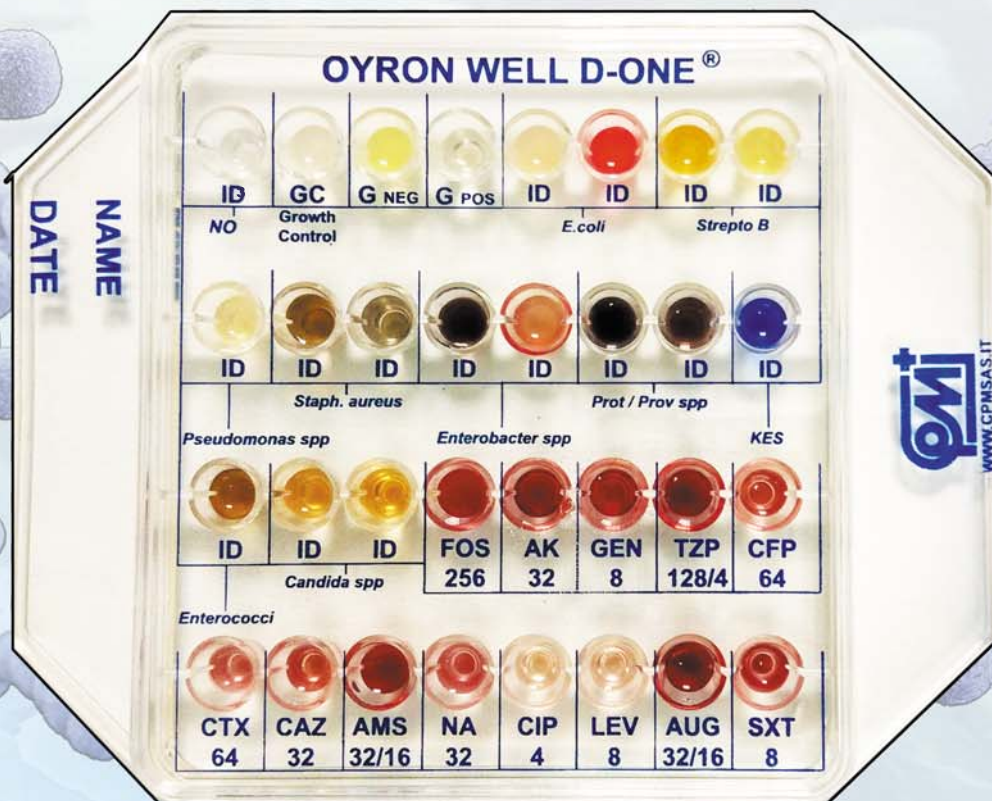

System for the presumptive identification and antimicrobial susceptibility test of:

Escherichia coli, *Streptococcus agalactiae*, *Pseudomonas spp.*,
Staphylococcus spp., *Staphylococcus aureus*, *Enterobacter spp.*,
Enterococcus spp., *Proteus/Providencia spp.*, *Pseudomonas spp.*,
KES group (Klebsiella, Enterobacter, Serratia), *Candida spp.*, *Candida albicans*.

The system also contains differential wells for Gram + and Gram - microorganisms and growth control well.

The tested sample is diluted in accordance with the procedure and inoculated directly into the wells of the plate.
Incubation is performed for 18-24 hours and readings are taken visually without additional equipment.

